### VITAL INJECTOR



#### SPECIFICATIONS

Needle Information	31 G, 9Pin Multi Needle Type
Needle Penetration Control	0~5 mm(Interval 0.2mm)
Syringe Type	1cc, 2cc, 3cc, 5cc
Input Power	AC 110~240V, 50/60Hz
Dimension (WxDxH mm)	Body : 388.8 X 231 X 133.4 Injector : 294.4 X 215 X 97.5
Weight	4.84Kg (Body + Injector)

**Eunsung** 



Mail: contacto@healthcorp.ec Celular: 0999230608 Innovative Technology For Your Medical Aesthetic Needs


# VITAL INJECTORS

The Next generation Vital injector3 will assist the skin to blossom into its full form

VITAL INJECTOR

## VITAL INJECTOR

Vital Injector3 is multi-use injector with a capacity of loading interchangeable 1cc~5cc syringes that is versatile to various medical use. Improved the structure of injectors and needles compared to the previous device, allowing hygienic treatment.

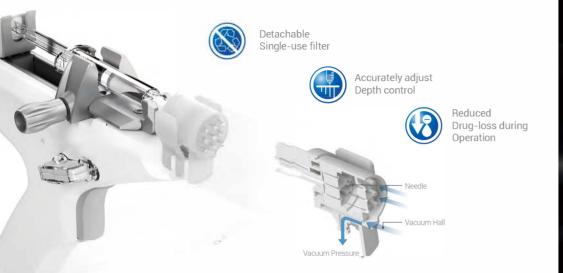



#### Disposable Multi Needle 9Pins in Round Tip

The circular designed Tip with 9 micro needles is reduced the pain sensation and treatment time. These Tips can only be used once, each pack is sterilized individually and allows hygienic and safe treatment.



## **Needle + Vacuum** penetration process





- 1 Controllable in range up to 5mm by 0,2mm interval.
- 2 Appropriate injection technique is dependent on the patient's skin and subcutaneous fat thickness
- 3 Drug is delivered at the target/suitable depth accurately.
- 4 Reduce treatment time by adapting the 31G, 9 Pin multi needle.

Mail: contacto@healthcorp.ec Celular: 0999230608

Minimum drug loss Maximum safety

2 A raised skin folds allows evenly insertion all of the 9 needles

and medication is not irritate or damage the beneath tissue.

1 The vacuum pressure lifts up skin fold




# An advanced anti-contamination system



#### Doble Filter System

An advanced anti-contamination system effectively protect patients and the device from contaminant. The closest filter with the injection site is designed as single-use that enhance hygienic of treatment. The secondary air filter has the Structure to easy disassemble.

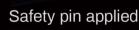

Designed to protect the user from accidental needlestick injuries.



Vital injector3 has a display screen that displays the injection information. User can set the Depth (needle insertion depth), Dose, Velocity and Vacumm level by simple operation of the UI.



- Syringe type setting (1, 2, 3, 5cc)
  - - Velocity setting (1~5 Level)
      Needle penetration Depth setting (0.0mm~5.0mm/0.2mm interval)
- Dose (Amount of injection) setting
- SHOT C/D (continuous/Dose mode)
- ✓ MODE setting SK (=skin) HR (=Hair, scalp)
- Convenient audible alarm signal
  - (Contents of syringe are almost consumed)